

Accurate Visualization During Laparoscopic Hysterectomy

apple MEDICAL

Colpo-Probe™ Vaginal Fornix Delineator



The Colpo-Probe includes an integral 10mm instrument guide with patented airtight seal ... added versatility for special procedures.



- **Reusable – easy to clean and sterilize, long useful life helps contain costs.**
- **Provides excellent anatomical visualization and accurate identification of structural boundaries.**
- **Aids in bladder retraction and lateral displacement of the ureters.**
- **Simple, straightforward technique helps reduce procedure time and improve efficiency.**
- **Maintains pneumoperitoneum during and after colpotomy and during vaginal vault closure.**


The Colpo-Probe is a simple, low cost instrument for laparoscopic hysterectomy that helps enhance operative precision while significantly reducing cost per procedure. Featuring a unique canted cup design, the device fits securely in the vaginal fornix, providing a clear, unambiguous backstop for the circumferential colpotomy. The extended barrel of the cup helps prevent loss of the pneumoperitoneum throughout the procedure.

Operating technique with the Colpo-Probe is simple, direct and easy to master. By applying pressure to the vaginal fornix, the cervical/vaginal junction is well presented allowing incision at the top of the vagina, as close to the cervix as possible. As a result, the instrument helps minimize "over-dissection" thereby optimizing maximal vaginal length and preserving uterosacral ligament support and nerve supply.

Best of all, the Colpo-Probe is both inexpensive and reusable. Fabricated of a medical-grade copolymer, the instrument is rugged, easy to clean and can be repeatedly sterilized by standard technique.

Accurate Visualization During Laparoscopic Hysterectomy

apple MEDICAL

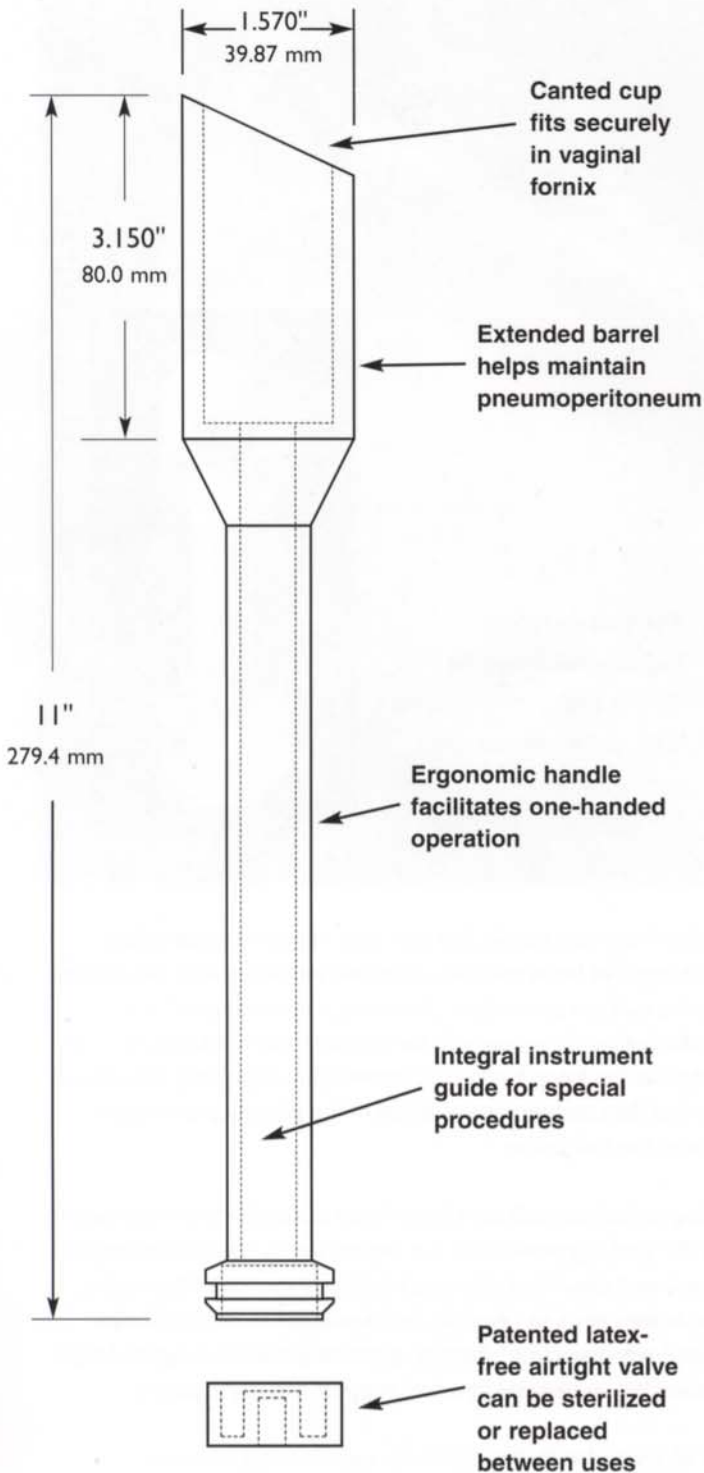

Colpo-Probe™ Vaginal Fornix Delineator

Ordering Information

Catalog No. 900-580

Colpo-Probe Vaginal Fornix Delineator

Includes one reusable instrument and 10 snap-on airtight valves.

Catalog No. 900-586

Replacement snap-on airtight valves.

pack of 10

Look to Apple Medical ...

Our products are designed by gynecologists specifically for gynecological procedures. Each item in our rapidly growing family of products is recognized for its ease of use, effectiveness and low cost.

This basic design concept enables us to respond to the needs of gynecologists for highly functional instruments that are simple and easy to handle during procedures. Apple Medical's innovative approach provides the physician with cost effective alternatives to the many expensive products currently offered to the gynecologist.

Look to Apple Medical Corporation as we continue to provide new instruments that meet the requirements of today's gynecologist and that embody the characteristics of our creative design philosophy: *Unique products ... simple solutions*

Call 800-255-2926 • Fax Orders 800-382-0349

Apple Medical Corporation

28 Lord Rd. Unit 135, Marlboro, MA 01752

Tel: 508-357-2700 or 800-255-2926 Fax: 800-382-0349

www.applemed.com

Unique Products. Simple Solutions.