



Lone Star Retractor System™ Gynecology Kit

...better than another pair of hands

For a wide range of
gynecologic surgical
procedures

Benefits

- Essential components for key procedures:
 - Urinary incontinence procedures including sling
 - Pelvic organ prolapse
 - Vaginal mesh revision/excision
 - Vesicovaginal fistula repair
 - Perineoplasty
 - Rectovaginal fistula repair
- Easy to use
- Better exposure
- Cost effective
- Easy to store, stock and pick






Lone Star Retractor System™ Gynecology Kit

At last, surgical access the way it should be. Unparalleled exposure. Under your control. This unique retractor system delivers uncluttered access using counter-traction with small hooks. The Result. A surgical site that you can easily adapt as the procedure progresses.

Ease and Precision: 3 Simple Steps

1. Place the Retractor Ring around the initial incision. Insert elastic stay hooks into tissue and retract in directly opposing quadrants. Slots on ring allow easy placement for equal and constant tension.
2. As dissection progresses, reposition hooks and add new ones as needed to achieve full, balanced retraction from any angle.
3. When closing, place stays along axis of the incision to create tension to facilitate precise staple or suture placement.

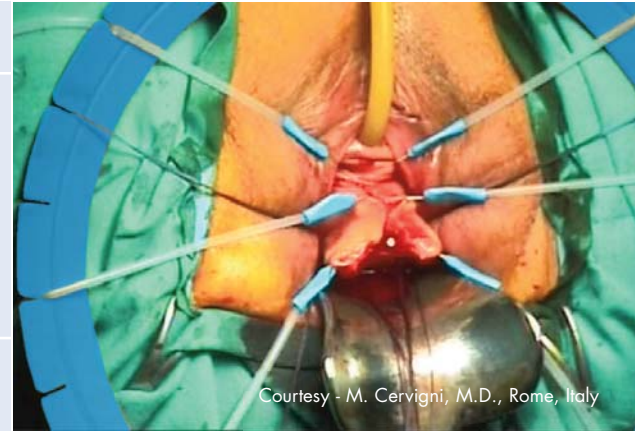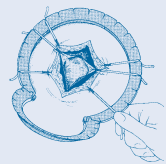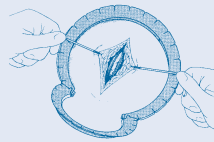

Courtesy - M. Cervigni, M.D., Rome, Italy

The Lone Star Retractor System facilitates mid-line incision of the anterior vaginal wall during cystocele repair.

Less Clutter. Better Exposure.
Effective Surgical Site
Management.

Get convenience without added cost!

Each kit contains one large snowman disposable retractor (31.8 cm x 18.3 cm – Model #3304G) and an eight-pack of 5 mm sharp hook elastic stays (model number 3311-8G) packed in a single-use sterile dispensing tray.

The cost is the same as if purchased separately, but the convenient packaging facilitates inventory management, making it easy to store, pick and present to the field.

Ordering Information:

Product Number	Description
3704	Each Lone Star Retractor System Gynecology Kit contains: <ul style="list-style-type: none">■ One large snowman disposable retractor ring, 31.8 cm x 18.3 cm – Model# 3304G■ One eight-pack of 5 mm sharp hook stays – Model# 3311-8G

Each kit is packaged in a single-use sterile dispensing tray. Sold in boxes containing 5 kits.

CooperSurgical

www.coopersurgical.com

Lone Star RETRACTOR SYSTEM™

95 Corporate Drive, Trumbull, CT 06611
800.243.2974 • 203.601.5200

Lone Star Retractor System™ is a trademark of CooperSurgical.

U.S. Patent Nos. 4,274,398; 4,430,991; 4,434,791; 5,769,783; 5,785,649; 5,964,697.

Form # 81451 • Rev. 09/08

© Copyright 2008 CooperSurgical, Inc