


Pipelle®

The most widely recognized name in endometrial sampling

Features and Benefits

- Consistently reliable
- Excellent patient acceptance
- Extraordinarily safe
- Rapid, simple technique
- Clinically-documented effectiveness
- Maximum sample yield
- Made in the U.S.A.

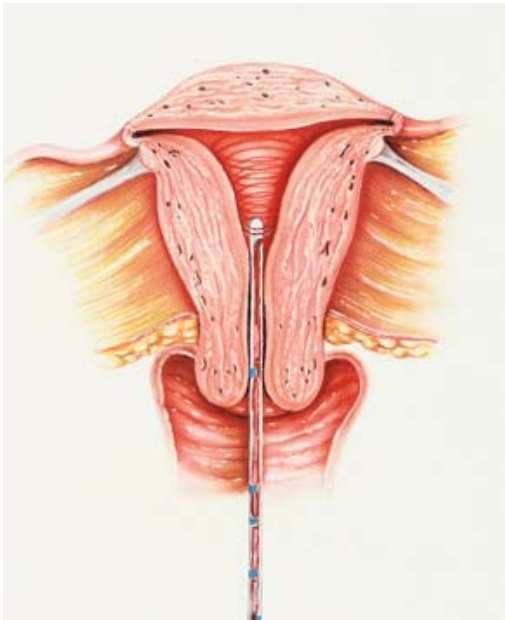

Pipelle[®]

Rapid, simple technique with proven clinical performance

- 3.1mm O.D. and flexibility allows for ease of insertion into most uteri
- Internal piston creates a vacuum for sample collection without any sharp or painful edges
- Procedure comfortably performed in as little as 30 seconds following basic patient preparation
- Completely self-contained with no noisy pump or clumsy syringes

Safe, Simple Sample Collection

Once within the uterine cavity, a vacuum is established within the Pipelle by manually retracting its internal piston. Sample collection is easily achieved by simply rotating the device constantly between your fingers as you slowly pass the Pipelle several times between the fundus and the internal os. The Pipelle provides unparalleled patient acceptance while achieving optimum uterine wall contact and maximum tissue yield.

Pipelle[®] is the “Gold Standard” in Endometrial Sampling. Utilizing a Simple Technique with Proven Clinical Performance and Reliability.

Created by a gynecologist for the gynecologist. Pipelle’s insightful design and exacting manufacture provide you with safe, simple, effective and cost-efficient performance that can be relied upon time, after time, after time.



Pipelle Ordering Information

Product Number: 8200

Packaging: Sterile, individually-packaged units of 25

CooperSurgical

CooperSurgical is the leading company dedicated to providing medical devices and procedure solutions that improve health care delivery to women regardless of clinical setting. Our company is fostering that position through expansion of its core businesses and introduction of advanced technology-based products which aid clinicians in the management and treatment of commonly seen conditions.

203.601.5200 | 800.243.2974 | coopersurgical.com | 95 Corporate Dr. | Trumbull, CT 06611
also at www.YouTube.com/CooperSurgical