

# Gynecath

## Uterine Infusion Catheter

Comfort and control  
in a single device

- **Comfortable**


- A soft, flexible, minimally invasive catheter minimizes discomfort during uterine exams
- The flexible plastic catheter, and latex-free balloon are gentle on delicate tissues

- **Easy to control**

- Singlehanded catheter control
- Once catheter is placed, there is no need for further manipulation
- Significantly reduces the need for a metal tenaculum

- **Easy to use**

- Sliding plastic sheath aids in catheter placement
- Latex-free balloon inflates easily once catheter is placed
- Balloon reduces fluid leakage, allowing optimal diagnostic viewing
- Available in a tray set to facilitate procedures

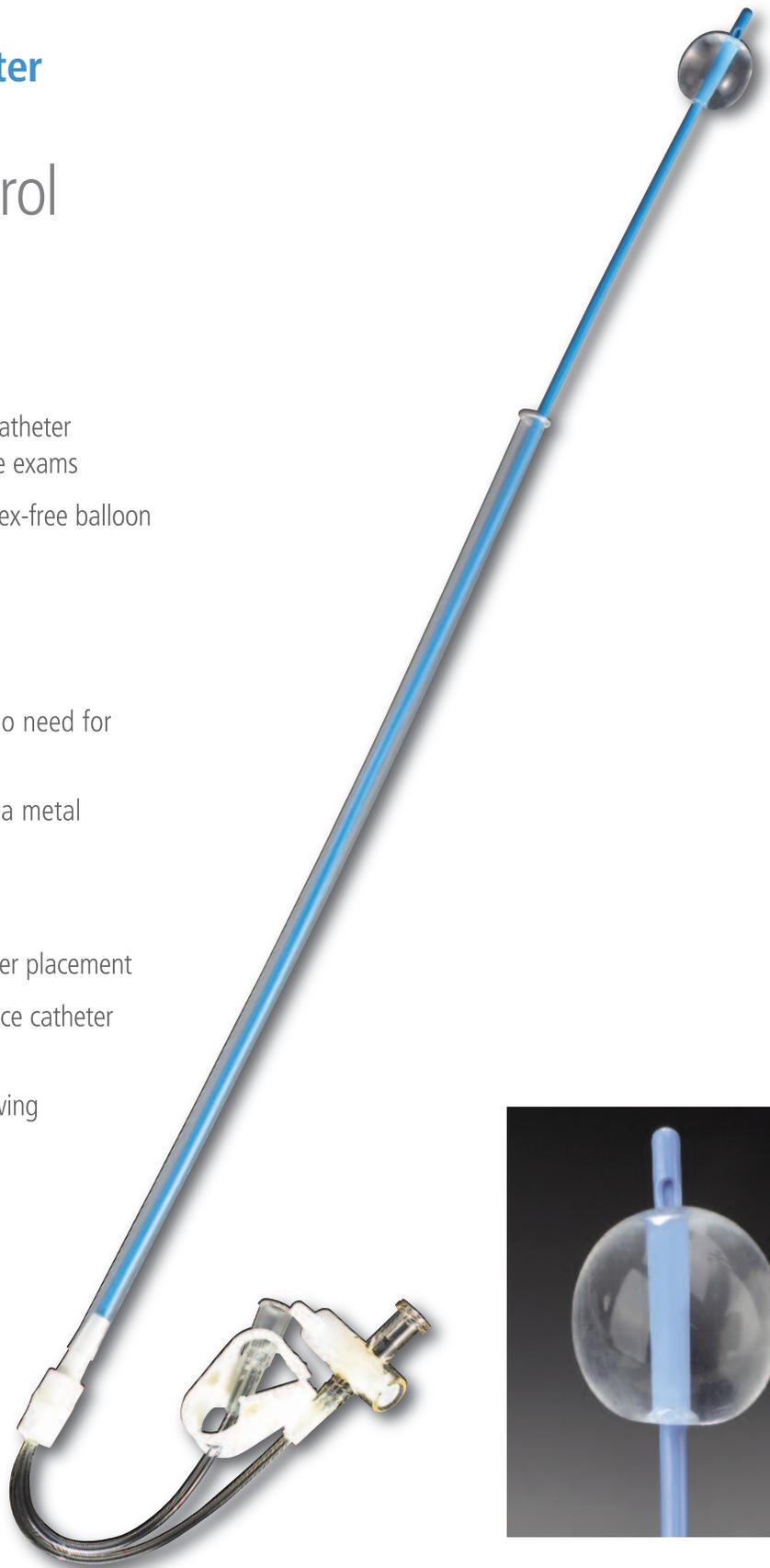

# Gynecath offers comfort and control in a single device

## Gynecath

The GyneCath Uterine Infusion Catheter is a minimally invasive, easy-to-use device that ensures a diagnostically enhanced, yet gentle, patient exam. Designed for use in saline-infusion hysterosonography and water-based hysterosalpingography.

## Enhanced Diagnostic Images

GyneCath is designed for use in saline-infusion transvaginal ultrasound and in radiologic exams using water-based contrast media. Once the balloon is inflated and fluid infused, the clinician's view of the uterine cavity and fallopian tubes is significantly enhanced.

## Product Specifications

Radiopaque 12-in. flexible 5 French and 7 French catheters have a closed, rounded distal tip with a side port proximal to the tip. The attached latex-free balloons have either a 3-cc or 6-cc capacity and are proximal to the side port. The dual lumen tube allows for balloon inflation and injection of the contrast media or saline solution through the side port.

## GyneCath Uterine Infusion Catheter

- Two sizes available:
  - 5 French for non-parous patients
  - 7 French for previously pregnant patients
- Available in tray sets for clinician convenience.



## Procedure Tray Sets

Each Procedure Tray contains:

- Fenestrated drape
- Disposable speculum
- Disposable cervical dilator
- Three swab sticks
- Sterile wrap
- Inner container for povidone-iodine solution
- One pack of povidone-iodine solution
- 36-in. extension tube
- 20-cc syringe with 18-gauge needle
- Lubricating jelly packet
- Sanitary napkin

## Gynecath Ordering information

Code	Product Description
Gynecath Catheters	
19610	Gynecath Catheter, 5 Fr, (10 per Box)
19611	Gynecath Catheter, 7 Fr, (10 per Box)
Gynecath Tray Sets	
19612	Gynecath Catheter, 5 Fr, (10 per Box)
19613	Gynecath Catheter, 7 Fr, (10 per Box)

## CooperSurgical

CooperSurgical is the leading company dedicated to providing medical devices and procedure solutions that improve health care delivery to women regardless of clinical setting. Our company is fostering that position through expansion of its core businesses and introduction of advanced technology-based products which aid clinicians in the management and treatment of commonly seen conditions.

203.601.5200 | 800.243.2974 | [coopersurgical.com](http://coopersurgical.com) | 95 Corporate Dr. | Trumbull, CT 06611  
Also on YouTube at [www.YouTube.com/CooperSurgical](http://www.YouTube.com/CooperSurgical)