



Your Trusted Source for Endometrial Sampling Devices

From the company dedicated to improving the lives of women, babies, and families.

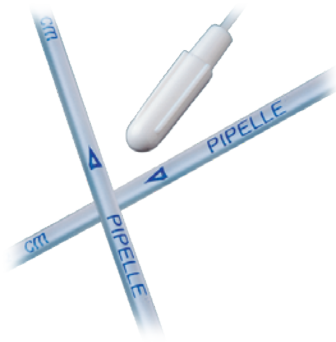

Pipelle® Endometrial Suction Curette

Single-use, sterile, disposable, suction curette for obtaining a histologic biopsy of the uterine mucosal lining or sample extraction of uterine menstrual content for microscopic examination or culturing. Clear, flexible, polypropylene sheath marked with colored, graduated markings from 4 cm to 10 cm, distance from the distal tip of the sheath to indicate the depth of insertion of the sheath into the uterus during use.¹

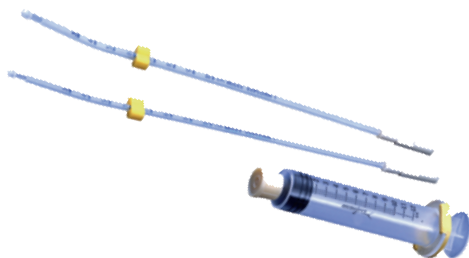

Uterine Explora™ Model I

Curette with Sharp Randall-Type Cutting Edge and Vacu-Lok Syringe

The MX120 (3 mm) and the MX124 (4 mm) curettes are used to obtain clearly differentiated endometrial tissue, usually without anesthesia, in an office setting. The flexibility of the nylon curette, the small diameter, and the sharp Randall-Type cutting edge allow better sampling from all areas of the endometrial cavity, including the cornual areas. The rounded end of this curette minimizes the possibility of perforation.³



Pipet Curet™ Endometrial Suction Curette

Single-use, sterile, disposable, suction curette for obtaining a histologic biopsy of the uterine mucosal lining or sample extraction of uterine menstrual content for microscopic examination or culturing.²



Uterine Explora™ Model II

Curette with Dual Hole Design and Vacu-Lok Syringe

The MX130 (3 mm) OD nylon cannula-curette is used to obtain endometrial tissue. The dual hole design of the MX130 curette obtains endometrial tissue via a scraping and peeling action. The dual hole design is particularly suitable for those patients who may have a low pain threshold. It is also the instrument of choice when a minimum of tissue is needed with the least disturbance to the endometrium. The ease of insertion and reduction in trauma enables the routine use of the Explora Curette as an outpatient or office procedure in the vast majority of cases, without anesthesia.⁴

Product Ordering Information

ENDOMETRIAL SAMPLERS

Pipelle® Endometrial Suction Curette

NUMBER	DESCRIPTION
8200	Pipelle Endometrial Suction Curette (25/Box)

Pipet Curet™ Endometrial Suction Curette

NUMBER	DESCRIPTION
MX140	3 mm (50/Box) – Rigid with greater integrity for stenotic Os
MX145	3 mm Soft Yellow (50/Box) – Pliable for molding catheter
MX150	4 mm (50/Box) – Multiparous patient

Uterine Explora™ Model I

Curette with Sharp Randall-Type Cutting Edge and Vacu-Lok Syringe

NUMBER	DESCRIPTION
MX120	Uterine Explora with Vacu-Lok Syringe, 3 mm, Model I (25/Box)
MX124	Uterine Explora with Vacu-Lok Syringe, 4 mm, Model I (25/Box)

Uterine Explora™ Model II

Curette with Dual Hole Design and Vacu-Lok Syringe

NUMBER	DESCRIPTION
MX130	Uterine Explora II with Vacu-Lok Syringe Model II (25/Box)

ADDITIONAL DISPOSABLES



Fundus Sound™ XM

NUMBER	DESCRIPTION
1196	Fundus Sound XM Disposable Uterine Sound (25/Box)



Os Finder™

NUMBER	DESCRIPTION
1176	Os Finder (50/Box)

To place an order, call us at: **800.243.2974**
or visit us online at: **www.coopersurgical.com**

Scan to Contact
a Representative



1. Pipelle® Endometrial Suction Curette IFU
2. Pipet Curet™ Endometrial Suction Curette IFU
3. Uterine Explora™ Model I IFU
4. Uterine Explora™ Model II IFU